




For rapid sterilization

EQUITRON®

TABLE TOP STERILIZER with drying

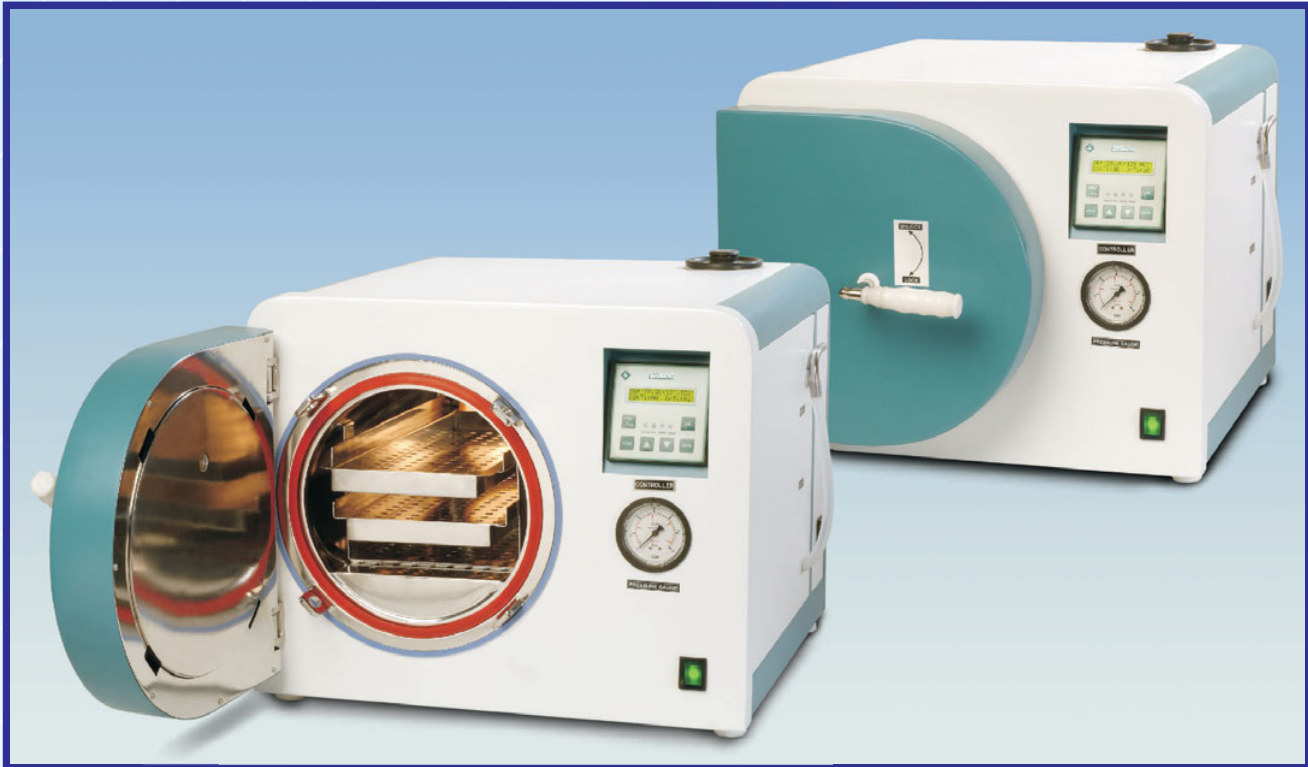

Applications:

- ❖ EQUITRON Table Top High Speed Sterilizers are ideal for all applications requiring rapid, reliable and total destruction of all living microorganisms.
- ❖ Ideal for Ophthalmic Surgeons, Dental Surgeons, Surgical Doctors, Polyclinics and Hospitals.

EQUITRON Advantages:

- ❖ Sterilizing, exhaust of water and steam drying are fully automatic.
- ❖ Requires a sterilizing time of just 7 minutes at 134°C (instead of 20 minutes at 121°C).
- ❖ A totally hassle and tension free sterilizer which can even be kept inside an operation theatre.
- ❖ Results reproducible because full cycle is automatic, requiring no human intervention.

Construction:

- ❖ Full stainless steel construction.
- ❖ Single lever pronged closure system with patented Silicone Gasket (registered design No. 208774.)

Standard Features:

- ❖ µP Controller with following parameters programmable by user: **Sterilizer Time:** up to 99 minutes. **Sterilizing Temperature:** up to 134°C. **Drying time:** up to 99 minutes.
- ❖ Pt100 Sensor for precise monitoring.
- ❖ End cycle buzzer and auto reset.
- ❖ Valves for water filling and exhaust.
- ❖ Industrial Grade energy efficient heater.
- ❖ Drying with door open facility.

Safety:

- ❖ Sensor open alarm and over temperature protection.
- ❖ Pressure interlock on door.
- ❖ Spring loaded safety valve for over pressure.

Model	Chamber Size	External Dimension W/D/H in mm	Tray	Heater	App. Shipping Details (Air) Net / Gr / Volume
#7471-16	Ø 230 x 400 mm (16 lit)	560 x 500+125 x 450	2 trays	1.2 kW	50 kg / 65 kg / 0.3 m ³

All electrical appliances work on 230 VAC, 50Hz, single phase.
Due to continuous improvement in design, the product supplied may have modified features.



MEDICA INSTRUMENT MFG. CO.

105, Vasan Udyog Bhavan, Sun Mill Compound, Lower Parel (West)
Mumbai 400013, INDIA • ☎ +91-22/6618-9999 • 📠 +91-22/6618-9996
✉ sales@medicainstrument.com • www.medicainstrument.com



0045_0907-1-1