



The sleek and heavy duty chiller

# equitron®

## CHILLED WATER CIRCULATOR - 25L

### Applications:

- ❖ Ideally suitable for chilling and condensing operations.
- ❖ Ideal companion for Blister Pack Machines, 20 lit Rotary Vacuum Evaporators, and for general production requirements of chilling in closed loop systems.

### EQUITRON Advantages:


- ❖ Environment friendly CFC free hermetically sealed **rotary** compressor.
- ❖ Sleek & compact design with easy maneuverability.

### Construction:

- ❖ Full stainless steel construction.
- ❖ Sturdy angular welded frame.
- ❖ Panels powder coated.
- ❖ Four heavy duty castors for easy movement.
- ❖ Outer size W/D/H (mm): 650 x 620 x 840.

### Standard Features:

- ❖ Suitable for closed loop systems.
- ❖ In-built Digital Temperature Controller cum Indicator.
- ❖ Temperature range upto +5°C (at ambient 25°C).
- ❖ Reservoir capacity of approximately 25 lit.
- ❖ Heavy duty water circulation pump of 180 W with 10 LPM flow rate at 15 met head.
- ❖ With low liquid level indication & automatic cut-off.

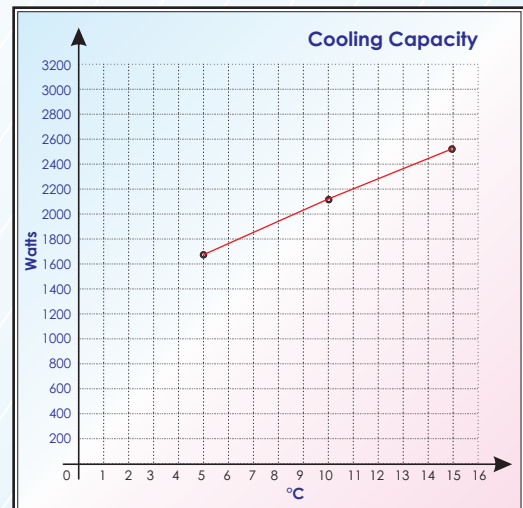

Model	Description	Cooling Capacity	Packing Dimensions (WxDxH) / Gr. Weight (WC)
#8502.025.000	Chilled Water Circulator - 25L	2500 W at 15°C	810 x 810 x 1090 mm / 129 kg

All electrical appliances work on 230V AC, 50Hz, single phase. Other Voltages available on request. Due to continuous improvement in design, the product supplied may have modified features.

**Our Other Products:** Autoclave Vertical: Standard / PAD / Fully Automatic, Autoclave Portable: Standard / Automatic, Autoclave SLL Series: Vacuum, Autoclave SLE series • Media Preparator • Table Top Autoclaves: Normal / Class B • Rotary Vacuum Evaporator • Vacuum Controller • Chilled Water Circulator • PTFE Diaphragm Vacuum Pump • BOD Incubator • Water Still • Water Bath: Stirred (Constant) / Serological / Concentric Rings / Refrigerated Circulating Bath • Round Bath: Clear Bottom / Immersion Series • Dry Bath • Colony Counter • Anaerobic Jar: SS / Polycarbonate • Pipette & Petri Dish Sterilizing Box • Test Tube Rack • ESR Westergren Apparatus • Slide Staining Stand / Rack / Trough • Slide Trays • Slide Cabinet • Test Tube Caps • Spreader • Cork Borer



**MEDICA INSTRUMENT MFG CO**

105, Vasan Udyog Bhavan, Sun Mill Compound, Lower Parel (West)  
Mumbai 400013, INDIA • ☎ +91-22/6618-9999 • 📠 +91-22/6618-9996  
✉ medica@medicainstrument.com • www.medicainstrument.com



0076\_1508-0-1